

SPLIT BULK STERILIZERS FOR EXISTING FACILITIES

REPLACING YOUR OUTDATED BULK AUTOCLAVE CAN BE EASY





Overcome your facility's egress limitations with Beta Star's exclusive sectional chamber designs. Chamber sections can be custom designed to fit through doorways, corridors and elevators, eliminating the need to disturb facility walls. All split bulks are custom designed, with section dimensions tailored around your egress requirements.

Sectional chambers are built as flanged pressure vessels in multiple sections based on egress and elevator limitations. Chambers are tested and certified to ASME Section VIII standards.

Once in your facility, the chambers are final assembled, welded and tested. Our piping skids are connected to the proper utilities and arrive with our exclusive EnviroVac® water conservation system.

Our contemporary HMI touch screen control provides users with simple steps to selecting, operating and monitoring sterilization cycles using a color, touch-sensitive, screen. Benefit from custom cycle programming, safety features, security settings and cycle monitoring.

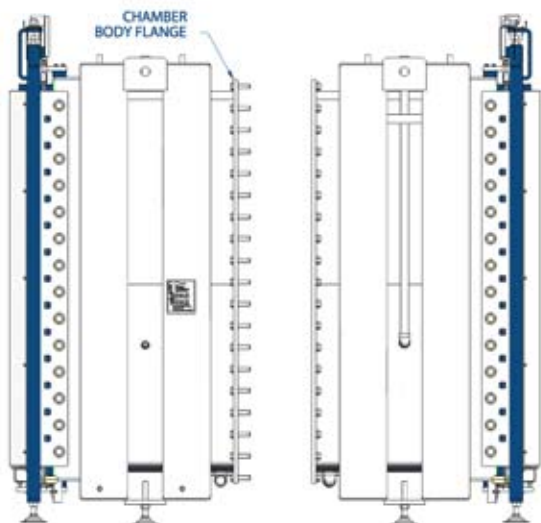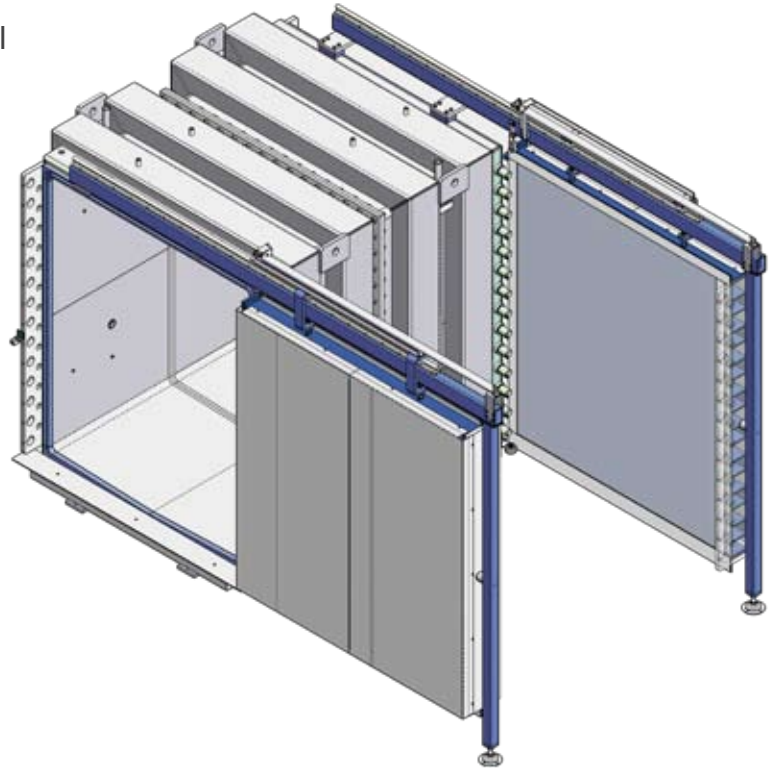

Eliminate interruptive and costly construction of your facility to creating an egress point. Vessel doors, door beams, and chamber sections are rigged into the facility separately through door openings and elevators.



SPLIT BULK FIT ASSEMBLED BEFORE SHIPMENT